



SABAMEDICARE

Point-of-care diagnostics · Lebanon

PRODUCT CATALOGUE

Multi-Function Monitoring Meter

Exclusive distributor in Lebanon

SabaMedicare SARL

For prices, demos and installation, contact us:

WhatsApp +961 71 490 154

info@sabamedicare.com

sabamedicare.com

Hazmieh, Baabda, Lebanon

Supplied in Lebanon by SabaMedicare SARL. Manufacturer specifications follow.



Multi-Functional Monitoring Meter



**7-in-1 meter,
for personalized health management**

Eaglenos Sciences, Inc.

Feature



7 in 1, multi-functional

Parameter	Sample Volume	Time
GLU 	0.8 μ L	5s
LAC 	0.8 μ L	10s
KET 	0.8 μ L	10s
UA 	1.5 μ L	10s
TC 	1.7 μ L	10s

HCT 
HB 

Wide application range

Parameter



Glucose



Lactate



Uric acid



β -Ketone



Total cholesterol



Hemoglobin



Hematocrit

Related diseases



Diabetes and complications




Hypertension



Hyperlipidemia



Renal injury

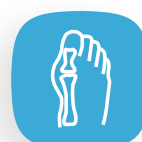


Hyperuricemia




Coronary heart disease

I Stay informed, stay healthy

Bluetooth connected and mobile APP

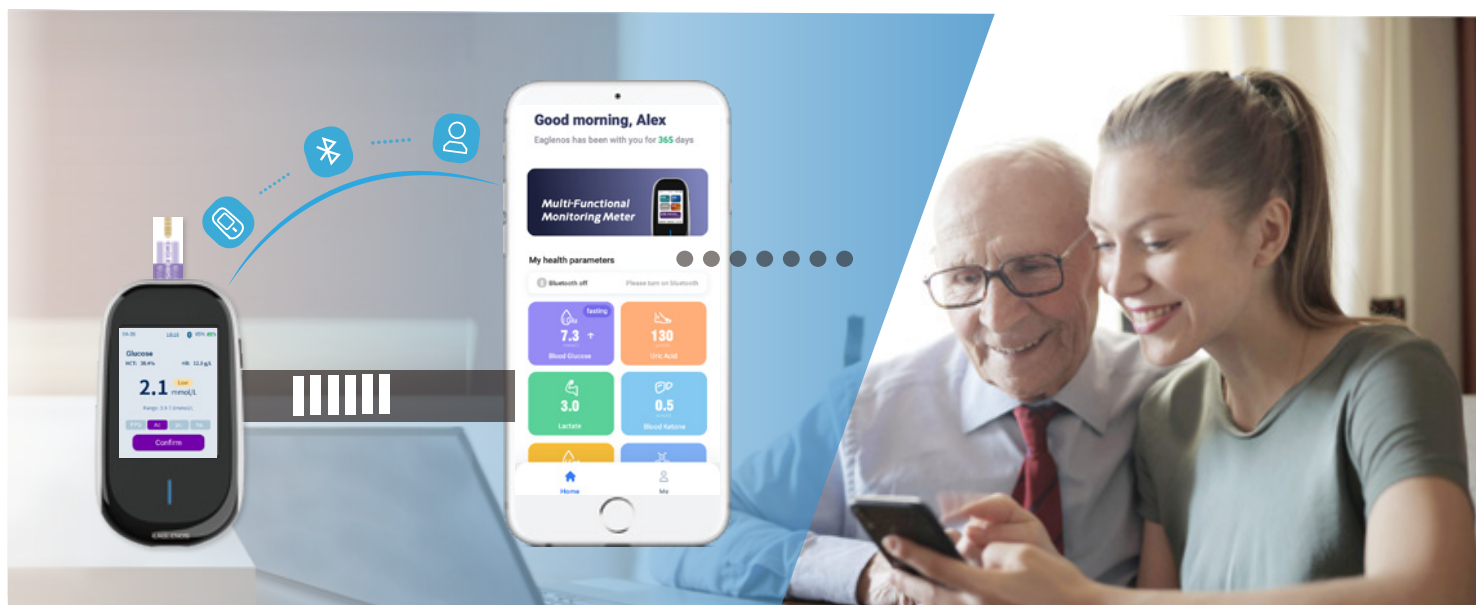

I 7-Gold-Electrode technology

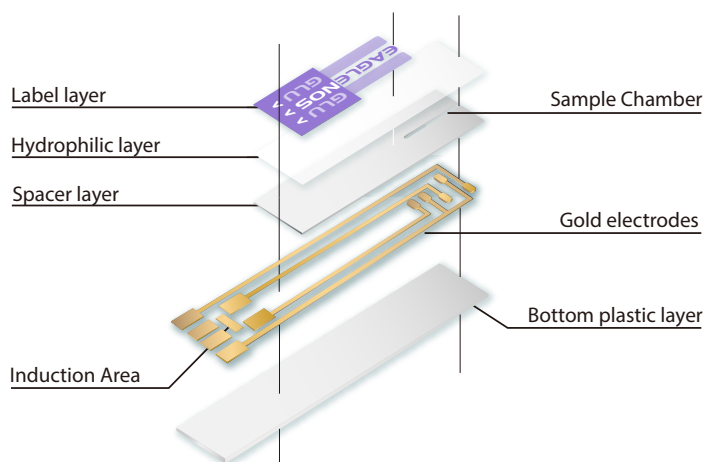

FAD-GDH technology to avoid blood oxygen interference



Enzymatic biosensor with high selectivity



Anti-interference enzyme chemistry design to prevent ascorbic acid interference



Sophisticated temperature compensation to ensure accuracy and consistency



Gold electrode, excellent sensitivity, conductivity and stability



Auto moisture-detection for damaged strip

I Ultra-wide HCT range



Blood glucose testing: HCT range 10% ~ 65%,
applicable for the vast majority of individuals

Adult males with Anemia

Non-pregnant adult women

Severe fetal anemia

Polycythemia (male)

Polycythemia (female)

Cardiovascular disease (female)

Cardiovascular disease (male)

Dialysis patients



HCT ranges of different populations

I Powerful partner for a healthy life



Smart identification

Auto-identification
of strip types



Smart coding

Smart coding using code chip with
calibration information



Large storage

Up to 7000 tests



User-friendly design

Built-in voice feature for
audible results



Position alert

Auto-alert when the strip
is in proper position



Convenience

Touch screen display



Eaglenos Sciences, Inc.

Tel: +86-400-019-0069

Website: www.eaglenos.com

Address: B2-2, Treehouse Headquarters Cluster, No.73 Tanmi Road,
Jiangbei New Area, Nanjing, Jiangsu 210044, China

*Please refer to the instruction manual for further information.



*For more information
please scan the code